CASADA
healthcare

< relaxmyway

QUATTROMED V

< relaxmyway

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< relaxmyway

Since the company was founded in the year 2000, the name CASADA has stood for high grade products in the area of wellness and fitness.

Today, CASADA – products are available worldwide. What distinguishes CASADA products is their continuous development, unmistakable design and their absolute focus on improving quality of life as well as their reliable and excellent quality.

Congratulations!

With the purchase of this massage device you have shown awareness of your health.

In order for you to enjoy the advantages of your massage device in the long term, we ask you to read and adhere to the instruction manual carefully.

We wish you a lot of fun with your personal Quattromed V.

Table of Contents

Safety Instructions.....	16
Scope of Delivery.....	17
Assembly / Startup	18
Troubleshooting	19
Specials.....	20
Remote Control.....	21
Storage / Maintenance	23
Specifications	23
Eu Declaration of Conformity.....	23
Warranty Conditions	24

Safety Instructions


Please read the following instructions carefully before putting your massage device into use to ensure trouble-free operation and optimum efficiency. Please retain these operating instructions for further use!

- The massage device conforms to recognised technical principles and the latest safety regulations.
- This device does not require any maintenance.
- Any possible repairs may only be carried out by authorised expert staff.
- Improper use and unauthorised repairs are not permitted for safety reasons and lead to loss of warranty.
- Never touch the power plug with wet hands.
- Please avoid contact of the device with water, high temperatures and direct sunlight.
- To avoid the danger of short-circuits, disconnect the power cord when not used for long periods of time.
- Do not use any damaged cables, plugs or loose sockets.
- If plugs or cords are damaged, they must be replaced by the manufacturer, a service representative or by qualified personnel.
- In case of malfunction, disconnect immediately from the mains.
- Misuse or incorrect use excludes any liability for damage.
- In order to avoid excessively stimulating the muscles and nerves, the recommended massage time of 30 minutes should not be exceeded.
- Each massage – even a hand massage – must be refrained from during pregnancy or if one or more of the following complaints are present in the massage area: recent injuries, thrombotic diseases, all kinds of inflammations and swellings, and cancer. It is recommended to consult a doctor before having a massage for the treatment of ailments and illnesses.
- Do not use this device while driving.
- If you are dependent on electrical aids e.g. pacemakers, please be sure to consult your doctor for medical advice before taking a massage.
- The packaging material supplied should not be used as a toy.

Scope of Delivery

- Power supply unit 24 volts
- Manual
- Seat section
- Backrest
- Remote control
- Allen key

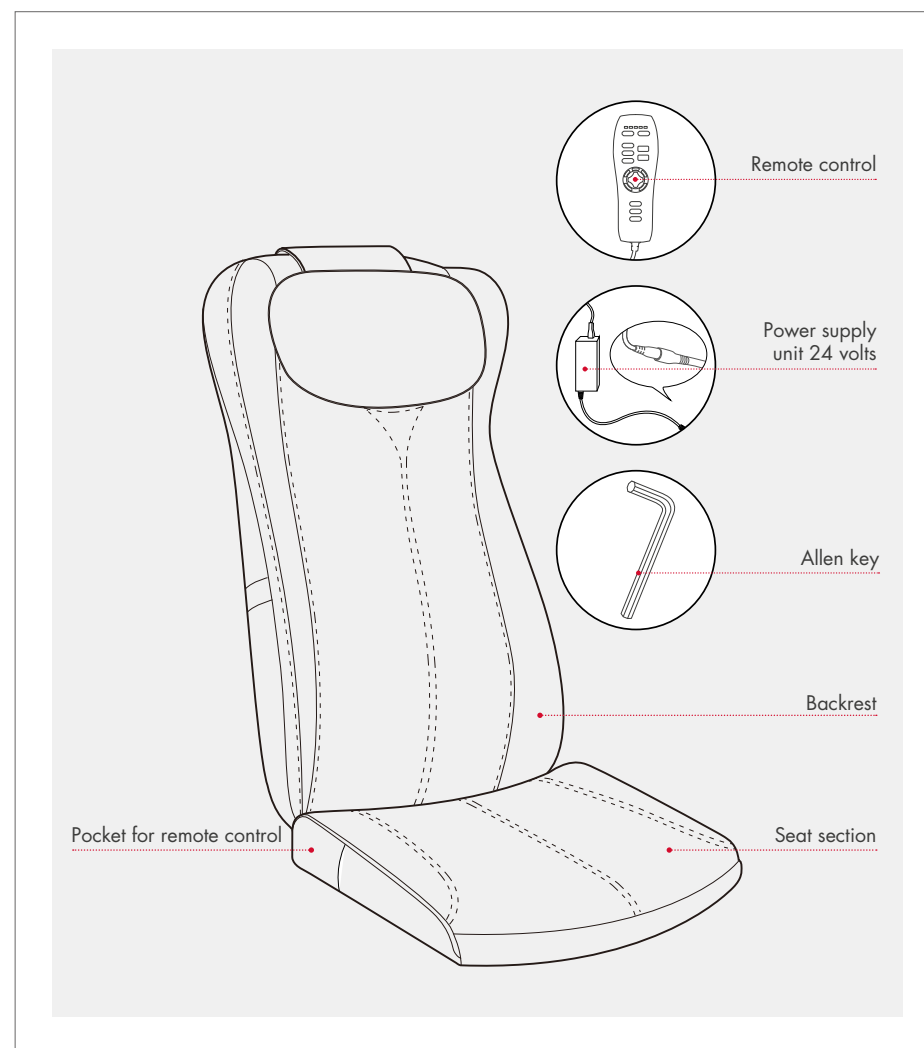
Assembly / Startup

1. Open the box, remove the power supply unit and place it within reach.
2. Take the black fabric bag out of the box carefully, open the zip fastener and remove the Quattromed V and the seat cushion packaged in protective film.
3. Remove the protective film.
4. The overlay/cover can be separated from the backrest via a Velcro fastener. Remove the protective cardboard between the two before first use. A padded backrest flap is provided for a gentler massage. Remove the backrest flap if you prefer more intense massage.
5. Back and seat section are not connected to each other at the time of delivery. For parallel use, remove the remote control from the pocket in the seat section and now attach both parts via the zip fastener.
(Illustration 1)
6. Then plug the cables together as shown. (Illustration 2)
7. Position the Quattromed V on a suitable chair or armchair. The massager is equipped with back strap that allow user to faster onto most chairs. Place the massager on a chair, wrap the two end of the strap at the back of the chair then secure with Velcro.
8. Finally, connect the power supply unit with the connection cable (Illustration 3), plug the power plug into the socket and turn the Quattromed V on using the red „ON/OFF“ button.


REMOVING LOCKING SCREW

Carefully loosen the screw with the attached Allen key on the back of the massage device. A screw is applied on the back of the product to protect the Quattromed apparatus during shipment. Remove the screw before use.


ATTENTION! The screw is not to be used again. Just remove it.

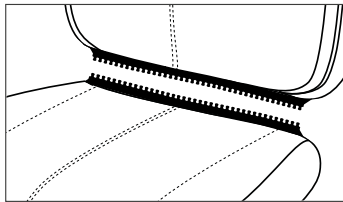(Illustration 1)

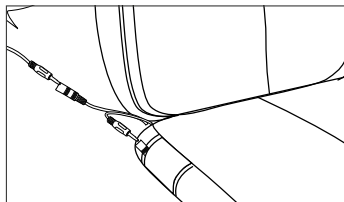(Illustration 2)

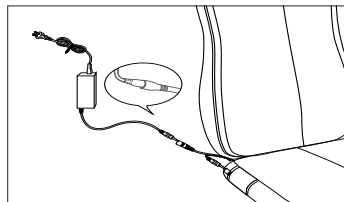(Illustration 3)

Troubleshooting


FAULT

- The device is connected to the mains, however neither the device nor the remote control operation work.
- The device is connected to the mains. The vibration in the seat is running and cannot be switched off. Remote control and massage do not work.
- The device is connected to the mains. Vibration massage does not work.


REMEDIAL ACTION

- Check whether all connections are correctly paired. Ensure that the power supply is fully connected.
- The device can overheat if used permanently. Therefore, turn it off, disconnect it from the mains power supply and let it cool down completely.
- Check whether the connectors are correctly paired for the vibration and power supply unit.


Specials

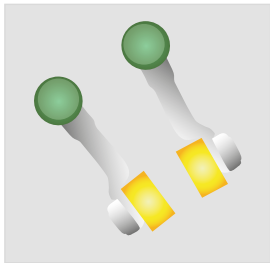The Quattromed V massages you with a total of four massage heads using the massage techniques „Shiatsu, Tapping, Kneading und Rolling“. Movable silicon massage heads as well as two fixed massage heads provide relief and mineral supply for the intervertebral discs. Massage heads, refined and enriched with certified jade care for a soothing and deep warmth massage.

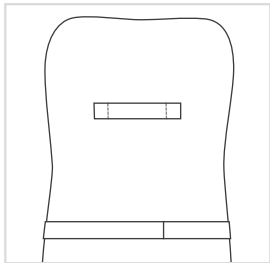The upper strap at the back serves as a support for arranging and positioning your Quattromed V and is not suitable for transportation needs. Please carry your massage device firmly with both of your hands while transporting. The lower Velcro strips serve as a fixation of the Quattromed V at the desired seating and offers a secured stand during the massage.

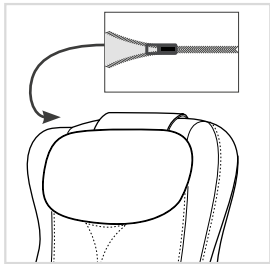The overlay/cover of your Quattromed V can be removed via the zipper and Velcro fastener on the backrest for cleaning or for a more intense massage.

For a vibration massage on your soles disconnect the seating part via the zipper and place it comfortably beneath your feet.

Specials

ROLLING

Pleasant alternation between relaxation and releasing the muscles – also perfectly suited for a final massage, for relaxation and recovery of the muscles.

KNEADING

Skin and muscles are either clasped and kneaded between the thumb and the index finger or with both hands. This massage technique is used in particular to ease tension.

VIBRATION

The Vibration setting moves parts of the body with equal vibrations, which penetrate deeply into the muscles. This stimulates the circulation in the muscles.

JADE HEATING

Jade is a natural store and conductor of heat. These high-quality jade stones, purposefully integrated in the massage heads dispense deep warmth and a relaxing effect.

AUTOMATIC

Specifically developed automatic programmes for the whole body, which offer a big range of different massage techniques. Additionally you can also turn on the heating function, if you'd like.

SPOT

The massage heads can massage a particular.

SHIATSU

Shiatsu (finger press) is a form of body therapy originating in Japan. The massage technique consists of soft, rhythmic, far-reaching stretches and rotations. The aim of Shiatsu is to simulate individual parts of the body and to mobilise the muscles.

WIDTH

The distance between the both massage heads can be set (narrow, medium and wide).

HEIGHT

Regulating the height limit of the massage technique in two steps.

3D MASSAGE

In this new technique, the massage robot not only moves upwards and downwards, but also forwards and backwards.

Remote control

- 1| **POWER**
Turn ON/OFF the power.
- 2| **HEIGHT**
When [Area]/[Auto] function is selected, press [Height] to select the lower limit of height. To exit the height limit function press [Height] again.
- 3| **AUTO PROGRAMME**
The massage device is equipped with 3 auto programmes. Press [Program], the unit is working during the cyclic switching of Mode 1 (Blue)/ Mode 2 (Green)/ Mode 3 (Red). The motor will preset first, and after finishing presetting, the unit starts the auto mode.
- 4| **HEATING**
The [Heating] function is activated at the beginning of following programmes [Tapping]/[Shiatsu]/[Auto]/[Area]. Press [Heating] button during massage to exit or switch on again the heating function.
- 5| **VIBRATION**
The [Vibration] can be switched on and off at any time during [Tapping] / [Shiatsu] / [Auto] / [Area] programme. Press [Vibration] button to start vibration function. The vibration will run in 1,5 min rhythm.
- 6| **SHIATSU**
Press [Shiatsu] button to start Shiatsu function. The default speed is middle. If the unit works in auto mode, it will exit from auto mode automatically. Press [Shiatsu] again to exit.
- 7| **TAPPING**
Press [Tapping] button to start tapping function. The default strength is middle. If the unit works in auto mode, it will exit from auto mode automatically. Press [Tapping] again to stop tapping function.
- 8| **SLOW – FAST**
The [Slow – Fast] function will only work when [Tapping] is selected. Press [Slow – Fast] button and the tapping strength will automatically cycle from weak to strong.
- 9| **FAST – SLOW**
The [Fast – Slow] function will only work when [Tapping] is selected. Press [Fast – Slow] button and the tapping strength will automatically cycle from strong to weak.
- 10| **WIDTH**
The [Width] function is only available with [Tapping] function. It does not work with [Shiatsu] function. Press [Width] button -1st -2nd -3rd step- the massage heads will enlarge from narrow to middle to wide and back. Press [Width] button -4th step- the device activates an continuously massage heads cycle from narrow to wide.
- 11| **FAST**
The tapping strength is classified into 3 levels: High-Middle-Low. The [Fast] button will only work when [Tapping] is selected. Press [Fast] button and the tapping strength will enhance to a higher level.
- 12| **SLOW**
The tapping strength is classified into 3 levels: High-Middle-Low. The [Slow] button will only work when [Tapping] is selected. Press [Slow] button and the tapping strength will decrease to a lower level.

Remote control

- 13| **LO**
The Shiatsu strength is classified into 3 levels: High-Middle-Low. The [LO Shiatsu Speed] will only work when [Shiatsu] is selected. Press [LO Shiatsu Speed] button and the shiatsu strength will decrease to a lower level.
- 14| **HI**
The Shiatsu strength is classified into 3 levels: High-Middle-Low. The [HI Shiatsu Speed] will only work when [Shiatsu] is selected. Press [Hi Shiatsu Speed] button and the shiatsu strength will enhance to a higher level.
- 15| **UP/DOWN**
Control keys of the massage heads for exact positioning on your back. During the localised massage, you can alter the massage area upwards or downwards by holding down this button.
- 16| **FULLBACK**
This automatic programme is for the entire back. Calming massage of the individual parts of the entire back with massage technique "Rolling". Optionally, the following functions can be activated [Tapping]/[Shiatsu]/[Heating]/[Vibration].
- 17| **UPPER**
This automatic programme is for the upper back area. Calming massage of the individual parts of the upper back with massage technique "Rolling". Optionally, the following functions can be activated [Tapping]/[Shiatsu]/[Heating]/[Vibration].
- 18| **LOWER**
This automatic programme is for the lower back area. Calming massage of the individual parts of the lower back with massage technique "Rolling". Optionally, the following functions can be activated [Tapping]/[Shiatsu]/[Heating]/[Vibration].

Storage / Maintenance

With regular maintenance you can increase the service life of the Quattromed V and ensure safe operation at all times.

1. Recommended storage: Dry, free of dust.
2. Keep liquids or open flames away from the Quattromed V.
3. Clean the device with an imitation leather care and with a dry soft cloth.
4. Please do not dip your Quattromed V into water or other liquids during cleaning.
5. Let the Quattromed V dry out completely after cleaning.
6. Ensure that no moisture penetrates the device via the buttons in the remote control.
7. Disconnect the Quattromed V from the main power supply when it is not in use.

Specifications Quattromed V

Backrest dimensions:	75 x 45 x 20 cm
Seat part dimensions:	7,5 x 45 x 42 cm
Weight:	8,5 kg
Voltage:	
Input:	110-240V~, 50/60 Hz
Output:	24V  , 2500mA
Nominal power:	
In operation:	60 w
Automatic runtime:	15 minutes

Certificates:

Eu Declaration of Conformity

Compliant with European regulatory standards EMC 2014/30/EU and IVD 2014/35/EU.

Warranty Conditions

The retailer guarantees a statutory period of warranty for the product described on the reverse. Warranty is valid from the date of purchase. The purchase date must be documented with proof of purchase. The retailer will repair or replace free-of-charge all items covered in Germany for which faults have been discovered. A prerequisite is that you send the faulty device with proof of purchase before the guarantee or warranty period expires. The warranty claim shall not be applicable if a fault is established which has arisen e.g. due to external influences or as a consequence of a repair or modification which was not carried out by the manufacturer or an authorised dealer.

The following points are not part of the warranty:

- Wear parts, e.g. abrasives, working synthetical parts, fabric covers, overlay / cover
- Engine sounds up to 72 db to 30 cm
- Wear caused by improper use
- Commercial use e.g. exhibition use, repeated use without intermediate breaks of 40 minutes
- Exposing the device to heat, humidity of liquids
- Damage due to sharp and/or sharp-edged objects
- Wanton destruction, misappropriation, overloading or accidental damage
- Damage and soiling by domestic animals
- Contamination caused by body or hair care products
- Damage caused by use of acids, solvent or alcohol agents

The warranty granted by the retailer is limited to repair and exchange of the product. Under this warranty, the manufacturer or the retailer shall have no further liability and are not responsible for damages which have arisen due to disregard of the operating instructions and/or improper use of the product.

Where upon examination of the product by the manufacturer it is found that the warranty request does not concern a defect covered by the warranty or that the warranty period has expired, the customer must bear the costs for the inspection and repair.

As manufacturer, Casada ensures the delivery of faultless goods to buyers/partners. Proven production faults are entitled to an exchange after registration and confirmation.

Legal note: All contents do not contain healing statements. The diagnosis and treatment of illnesses and other physical disorders require a treatment by a physician, non-medical practitioner or therapist. The statements are strictly informative and may not be used as a substitute for a professional medical treatment. Every user is beholden to ascertain by accurate examination and, if required, after consultation of a specialist, if the application is the right one for their specific cases. Each dose or application takes place at the risk of the user. We dissociate ourselves from any healing statements or promises.

Casada International GmbH
Obermeiers Feld 3
33104 Paderborn, Germany
info@casada.com
www.casada.com

Copyright pictures and text are subject to the copyright of CASADA International GmbH and may not be used without explicit confirmation.

Copyright © 2017 CASADA International GmbH.
All rights reserved.

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